

ALDOLASE - CONTROL SET

Animal Control sera for Aldolase determination in serum and plasma

INTENDED USE

Concentrations and activities of components are studied for an excellent use even in automatic system of clinical chemistry, even in manual procedures.

PRESERVATION AND STABILITY

Preserve the product at 2-8°C, stable until expiration date on the package.

Once opened controls are stable for 5 days at 2-8°C.

MATERIALS REQUIRED BUT NOT SUPPLIED

Automatic Micropipette, normal laboratory equipment, distilled or deionized water.

REAGENTS PREPARATIONS

Open carefully one vial of liophile, trying to avoid loss of material, and with the enclosed pipette add exactly 1.0 ml of deionized water, for this is better employing an automatic pipette.

Close with attention and let the liophile recovers for 30 minutes, shake carefully trying to avoid foam formation. Bring controls at room temperature before use.

PRECAUTIONS

- Product contains Sodium Azide 0.095% as preservative.
- Avoid contact with skin and mucous membrane.
- Use the normal precautions for lab behaviour.
- For in vitro use only.

WASTE DISPOSAL

Product is intended for professional laboratories. Waste products must be handled as per relevant security cards and local regulations.

PACKAGING


CODE CC03330

Control Level 1 1 x 1 ml (liophile)
Control Level 2 1 x 1 ml (liophile)

REFERENCES

Archives data.

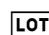

CALIBRATORS VALUES

 B030.26*

 2028/01

COMPONENT / METHOD	T (°C)	VALUE	RANGE	UNIT
--------------------	--------	-------	-------	------









Kinetics UV Aldolase determination in serum and plasma

CONTROL LEVEL 1	T (°C)	VALUE	RANGE	UNIT
 T013.26A 2028/01	37°C	10.06	8.05 – 12.08	U/l
CONTROL LEVEL 2	T (°C)	VALUE	RANGE	UNIT
 U013.26A 2028/01	37°C	20.82	16.65 – 24.98	U/l

MANUFACTURER

LTA s.r.l.
Via Milano 15/F
20041 Bussero (Milan) ITALY
Tel: +39 02 95409034
Fax: +39 02 95334185
e-mail: info@LTAonline.it
Website: <http://www.LTAonline.it>

SYMBOLS

-  Only for IVD use
-  Lot of manufacturing
-  Code number
-  Storage temperature interval
-  Expiration date (year – month)
-  Warning, read enclosed documents
-  Read the directions
-  Biological risk

Mod. 01.06 (ver. 1.1 - 17/07/2006)

